Embrace Prevention...

Giving the tooth what it needs

Fills superficial, non-carious enamel lesions (white spots)



CXP[™] = **X**ylitol-coated **C**alcium and **P**hosphate for unsurpassed fluoride release

Releases 10 times more fluoride than the leading brand.

There's more to Embrace than fluoride release


- Pleasing taste ensures patient compliance
- Contains bioavailable calcium, phosphate and fluoride
- Does not separate no mixing required
- Ensures predictable, uniform dose

Code	Description
FV50	50 x 0.4mL with brush
FV200	200 x 0.4mL with brush
FVT	Tube, 12mL
FVX100	100 x 0.4mL (no brush)





By incorporating xylitol-coated calcium and phosphate in a permeable resin matrix that does not separate, Pulpdent has developed a sustained time-release varnish with uniform dosage that delivers **10 times more fluoride than the leading varnish brand.** The pleasing taste ensures patient compliance.

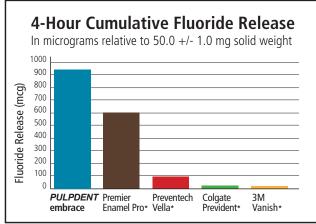

*Not a trademark of Pulpdent Corporation

THE DENIAL ADVISOR

Yapp R, Powers JM. Fluoride Ion Release from Several Fluoride Varnishes. Dent Advis Res Rpt 45:1, March 2012.

Bioavailable Fluoride, Calcium and Phosphate with Xylitol

Not only does Embrace Varnish release more fluoride in 4 hours, it also releases bioavailable calcium and phosphate ions, the essential building blocks of teeth. The xylitol coating prevents the calcium and phosphate salts from reacting until they come in contact with saliva. Saliva dissolves the xylitol and releases the calcium and phosphate ions, which react continuously in saliva with the fluoride ions to form protective fluorapatite on the teeth.

Thin-to-Win Application of Embrace Varnish




Dry the teeth with a gauze pad. It is best if the teeth are clean, but prophylaxis is not necessary.



Pinch varnish foil pack to push contents back from tear line, and tear at slit to open varnish pack.



Dispense varnish on glove or a pad. Pick up varnish with brush. (Embrace does not separate and requires no mixing.)

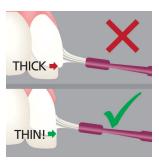

Apply a THIN coat of Embrace Varnish with one horizontal swipe of the brush. Apply Embrace Varnish *Thin-to-Win*.